




**Rail Mounting Kit with triple hook, factory attached, powder coated metal wire hook with integrated mounting plate and adapter**

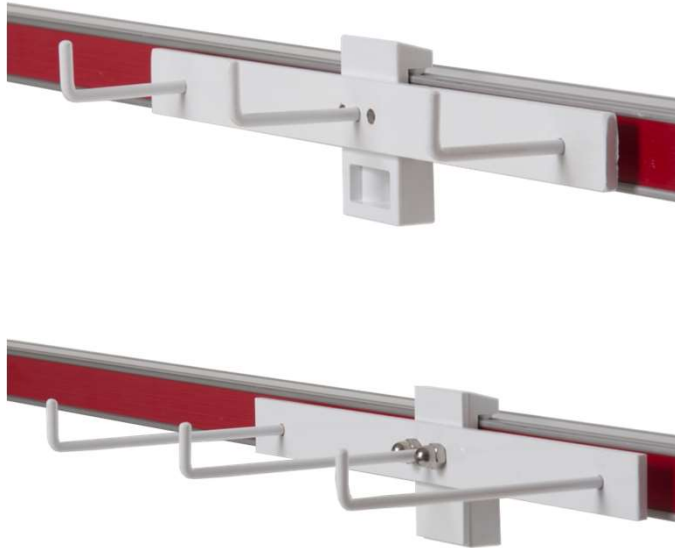

**Item Numbers:**

**SM302-0.75-PM (3/4" hook)**

**SM302-1.5- PM (1.5" hook)**

**SM302-2-PM (2" hook)**

**SM302-4-PM (4" hook)**

**JSN Number: not established**

**Dimensions:**

**12"W x 2.5"H**

**x length of hook"D, Weight: 2lbs**

**Item Description:**

- Mounts to "EVOLUTION," Fairfield-style Medical Equipment Management Rail using a snap action adapter with positive snap attachment lock feature.
- Integrated Allen screws provide the ability to lock the adapter in place on rail.
- Standard adapters are plastic with embedded anti-microbial compound.
- Adapter Plates are aluminum with an anti-microbial powder coat painted finish (AMF).
- Some accessories which require heavy-duty capacity will be provided with extruded aluminum adapter(s) and/or wall stand offs for added strength.
- Some plates have multiple holes to attach future products from different vendors without the need to purchase new plates each time.
- Some plates allowing adaptability to user demand, new codes or new vendor products without going to the wall and adjusting the rail itself.
- Alternate mounting kits available for other manufacturer's devices. Please contact your Sales Representative.
- Accessories now available for Modular, Hill-Rom/Wittrock and AMICO Style Rails.